

## Knotless FiberTak® Biceps Implant System

The Knotless FiberTak Biceps Implant System offers a knotless, all-suture anchor solution to perform an onlay proximal biceps tenodesis at any location, depending on surgeon preference. The implant uniquely integrates two individual tensionable knotless mechanisms within a single FiberTak sheath. The blue FiberWire® repair suture passes around the biceps tendon, while the white suture loop is retrieved through the tendon. The blue repair suture is retrieved through the loop and both are tensioned to complete the tenodesis.




Knotless FiberTak Biceps Implant System

### PRODUCT OVERVIEW

Each implant system contains a Knotless FiberTak Biceps anchor, a drill guide, a 3 mm drill pin, and a Loop 'N' Tack™ SwiftStitch™ suture passer.


### KEY FEATURES AND BENEFITS

- Two individual tensionable knotless mechanisms to secure the biceps tendon without whipstitching
- Versatile onlay option for any proximal biceps tenodesis location

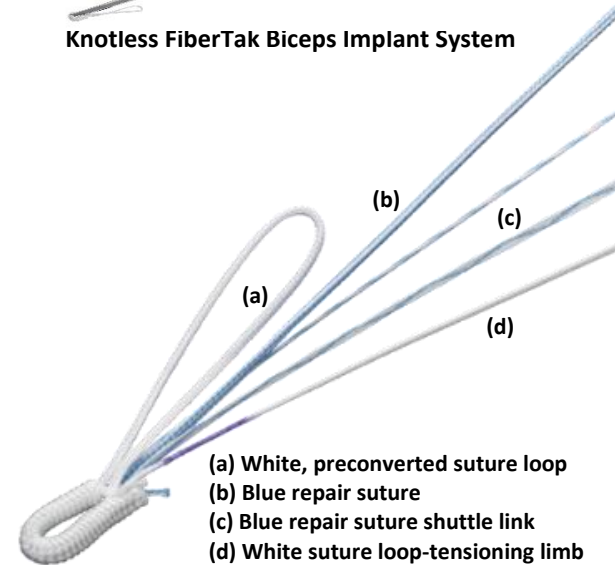

(a) White, preconverted suture loop  
(b) Blue repair suture  
(c) Blue repair suture shuttle link  
(d) White suture loop-tensioning limb


### SALES TOOLS AND RESOURCES

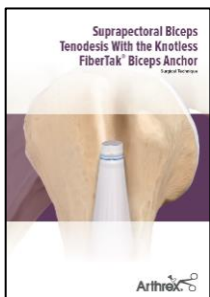


LT1-000251-en-US



AN1-000450-en-US

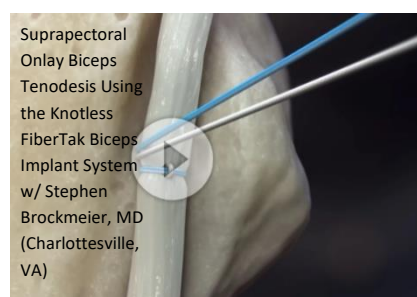

VID1-003955-en-US

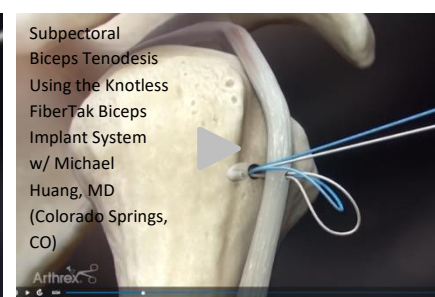

VID1-004078-en-US

### ORDERING INFORMATION

Part Number	Description	Stocking	LSP
<b>AR-3688</b>	<b>Knotless FiberTak Biceps Implant System</b> Knotless FiberTak Biceps Implant 3 mm Drill Drill Guide w/ Sharp Obturator Loop 'N' Tack SwiftStitch Suture Passer	<b>\$990</b>	<b>\$1650</b>
<b>AR-3687</b>	Knotless FiberTak Biceps Implant	<b>\$575</b>	<b>\$959</b>

Allen Holowecy

Sr. Product Manager, Shoulder & Elbow