


## 3.9 mm SwiveLock® Anchor


The 3.9 mm SwiveLock anchor leverages the proven thread and larger eyelet design of our most popular SwiveLock anchors but in a smaller footprint.

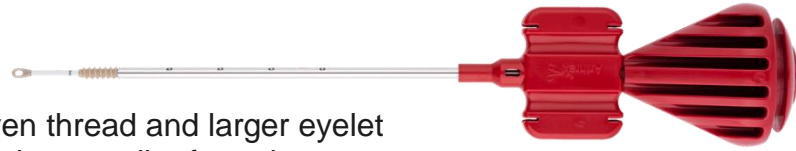

### PRODUCT OVERVIEW

The 3.9 mm SwiveLock anchor's **large eyelet size** makes it easy to use with FiberTape® suture, and its **deep threads** allow for strong fixation.

This combination makes it a great option where pathology requires additional fixation points but where space is limited (ie: Knotless SpeedBridge™ medial row, etc), or in small spaces where a deeper thread for robust fixation with tape suture is preferred (ie, *InternaBrace*™ applications, Loop 'N' Tack™ biceps repair, etc).

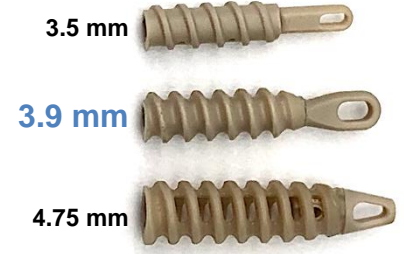

### KEY FEATURES & BENEFITS

- **Small size** combined with **large suture eyelet** makes it ideal for tight spaces
- Large suture eyelet **easily accommodates two FiberTape or three SutureTape sutures**, and auto releases from the driver as the anchor body advances
- **Hexalobe drive** (new design) for high insertional torque in hard bone
- **Instruments:** spade tip drill (packaged sterile, limited reusable), punch, and tap (for hard bone)



Hexalobe cross section




Spade tip drill (3.2 mm)



Punch (3.2 mm)



Tap (for hard bone)

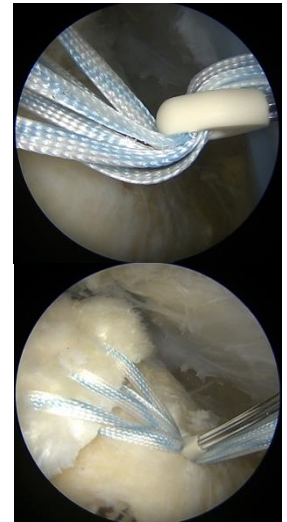

### SALES TOOLS & RESOURCES

**SwiveLock iBook:** excellent sales resource for scrub sink presentations and office visits  
[UE1-0456-EN](#)

A biocomposite version with vents is in development and will be available soon. Watch for procedure-specific kits coming in the future (Achilles SpeedBridge repair, etc).

### ORDERING INFORMATION

Part Number	Description	Stocking	LSP
<b>AR-2326PSLC</b>	PEEK SwiveLock Anchor 3.9 × 17.9 mm	<b>\$375</b>	<b>\$536</b>
<b>AR-2326BCC</b>	BioComposite SwiveLock Anchor 3.9 × 17.9 mm	<b>ETA FY20</b>	<b>QTR3</b>
<b>AR-2326P</b>	Punch for 3.9 mm SwiveLock Anchor	<b>\$190</b>	<b>\$272</b>
<b>AR-2326PTB</b>	Tap for 3.9 mm SwiveLock Anchor	<b>\$495</b>	<b>\$708</b>
<b>AR-2326D</b>	Spade Tip Drill for 3.9 mm SwiveLock Anchor	<b>\$105</b>	<b>\$150</b>

Loren Crook

Product Manager, Shoulder & Elbow

